

LENNTECH

info@lennotech.com Tel. +31-152-610-900
www.lennotech.com Fax. +31-152-616-289

Zeta Plus™ 1MDS Series Filter Discs and Cartridges

Zeta Plus 1MDS filter media is designated by the United States Environmental Protection Agency (Information Collection Rule, EPA/600/R-95/178, April 1996) to determine the level of total culturable virus in drinking water.

Zeta Plus 1MDS filter media discs and cartridges provide a convenient and cost-effective method for capturing and concentrating waterborne enteric viruses. Equally useful in the laboratory and in the field, Zeta Plus 1MDS products are suitable for sampling large water sources such as lakes, rivers, pools and sanitation effluent.

Zeta Plus 1MDS media uses electrokinetic adsorption to capture viruses. In this system, the unique electropositive charge of the media allows sample testing at the prevalent pH, eliminating all post-collection sample preparation such as chemical modification. Zeta Plus 1MDS also greatly reduces precipitation and related filter clogging caused by salts added to water sources.

Experiments involving adsorption rate, sample volume and virus viability have demonstrated the efficacy of Zeta Plus 1MDS products. Trials show that the medium is capable of processing large samples at high flow rates and that adsorbed virus remains viable up to 300 hours, eliminating the need for immediate sample processing.

Performance

Virus Elution

Adsorbed viruses can be recovered effectively from Zeta Plus 1MDS media using an alkaline buffered protein solution, such as a 1.5% beef extract in an aqueous glycerine solution at pH 9.5.

Stability

It is important that the concentration media have little or no effect on the virus viability, because it is not always convenient to elute collected virus in the field. When Zeta Plus 1MDS discs are stored at 4 °C, polio virus adsorbed to the media retain their infectious nature for up to 300 hours (12.5 days) with little appreciable loss.

Features and Benefits

Efficient capture & concentration of waterborne viruses

- Allows sample testing at prevalent pH

Adsorbed virus remains viable for up to 300 hours

- Eliminates need for immediate sample preparation

Media available in cartridges and discs

- Cost-effective and convenient use in the laboratory and in the field



Applications

Isolation of viruses in samples from:

Drinking Water

- Surface or ground supplies

Recreational Water

- Pools, lakes and rivers

Waste Water

- Sewage or industrial waste, including reclaimed water for crop irrigation

The 3M logo, consisting of the letters '3M' in a bold, red, sans-serif font.


Zeta Plus™ 1MDS Ordering Guide

1. Cartridges

Basic Cartridge Number	Nominal Cartridge Length	Media Formulation
45144	01 = 9 ¾ inch (24.8 cm) Nitrile Gasket 02 = 10 inch (25.4 cm) Nitrile Gasket	1MDS

2. Discs

Basic Disc Number	Disc Size	Media Formulation
64085	01 = 47 mm 02 = 90 mm 03 = 142 mm 04 = 293 mm	1MDS

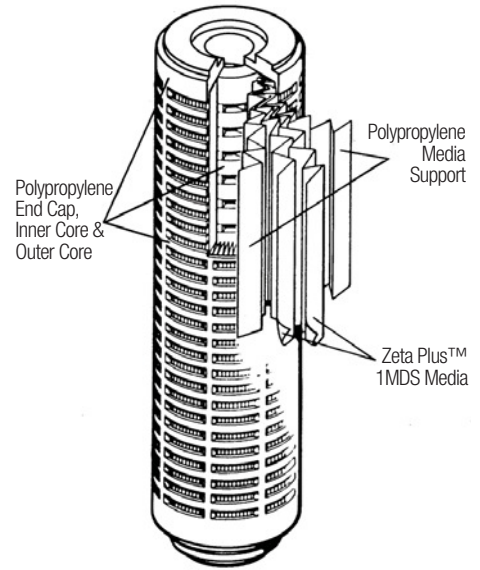

Construction

Zeta Plus 1MDS is an electropositive, charge-modified, glass and cellulose medium. Cartridges are constructed using two layers of pleated medium and are available in 10 inch (25.4 cm) and 9.75 inch (24.8 cm) double open-ended (flat nitrile gasket) configurations. Each cartridge contains approximately 0.39 m² of media and the cartridge may be sterilised using chemical or heat methods (autoclave to 121 °C or 1 hour). When using Zeta Plus 1MDS discs, it is recommended that two layers be used for optimal adsorption. As physical properties of different virus types and environmental conditions vary, it is recommended that the optimum system for the use of 1MDS be established for each laboratory. For more information, contact your local 3M Purification distributor.

Important Notice

3M Purification MAKES NO WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. Since a variety of factors can affect the use and performance of a 3M Purification product in a particular application, some of which are uniquely within the user's knowledge and control, user is responsible for determining whether or not the 3M Purification product is fit for a particular purpose and suitable for user's method of application.

Limited Warranty

3M Purification warrants it this product to be free from defects in material and workmanship during normal use for a period of one (1) year from the date of shipment from the factory. If the Product(s) is (are) defective within this warranty period, your exclusive remedy and 3M Purification's sole obligations shall be, at 3M Purification's option, to replace or repair the Product(s) or refund the original purchase price of the Product(s). This warranty does not apply to failures that result from abuse, misuse, alteration or damage not caused by 3M Purification or failure to properly follow installation and use instructions.

Limitation of Liability: 3M Purification will not be liable for any loss or damage arising from the use of the Product(s), whether direct, indirect, special, incidental, or consequential, regardless of the legal theory asserted, including warranty, contract, negligence or strict liability. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights and you may have other rights which vary from state to state, or country to country.

Zeta Plus and 3M are trademarks of the 3M Company used under license.



LENNTECH

info@lenntech.com Tel. +31-152-610-900
www.lenntech.com Fax. +31-152-616-289

Data may be subject to change without further notice.

v.2 - DOC10180
© 3M 2010. All rights reserved. LITCTZP1MDS1.E - 1110

Zeta Plus™ 1MDS Series

Filter Discs and Cartridges

Zeta Plus 1MDS filter media is designated by the United States Environmental Protection Agency (Information Collection Rule, EPA/600/R-95/178, April 1996) to determine the level of total culturable virus in drinking water.

Zeta Plus 1MDS filter media discs and cartridges provide a convenient and cost-effective method for capturing and concentrating waterborne enteric viruses. Equally useful in the laboratory and in the field, Zeta Plus 1MDS products are suitable for sampling large water sources such as lakes, rivers, pools and sanitation effluent.

Zeta Plus 1MDS media uses electrokinetic adsorption to capture viruses. In this system, the unique electropositive charge of the media allows sample testing at the prevalent pH, eliminating all post-collection sample preparation such as chemical modification. Zeta Plus 1MDS also greatly reduces precipitation and related filter clogging caused by salts added to water sources.

Experiments involving adsorption rate, sample volume and virus viability have demonstrated the efficacy of Zeta Plus 1MDS products. Trials show that the medium is capable of processing large samples at high flow rates and that adsorbed virus remains viable up to 300 hours, eliminating the need for immediate sample processing.

Performance

Virus Elution

Adsorbed viruses can be recovered effectively from Zeta Plus 1MDS media using an alkaline buffered protein solution, such as a 1.5% beef extract in an aqueous glycerine solution at pH 9.5.

Stability

It is important that the concentration media have little or no effect on the virus viability, because it is not always convenient to elute collected virus in the field. When Zeta Plus 1MDS discs are stored at 4 °C, polio virus adsorbed to the media retain their infectious nature for up to 300 hours (12.5 days) with little appreciable loss.

Features and Benefits

Efficient capture & concentration of waterborne viruses

- Allows sample testing at prevalent pH

Adsorbed virus remains viable for up to 300 hours

- Eliminates need for immediate sample preparation

Media available in cartridges and discs

- Cost-effective and convenient use in the laboratory and in the field



Applications

Isolation of viruses in samples from:

Drinking Water

- Surface or ground supplies

Recreational Water

- Pools, lakes and rivers

Waste Water

- Sewage or industrial waste, including reclaimed water for crop irrigation


Zeta Plus™ 1MDS Ordering Guide

1. Cartridges

Basic Cartridge Number	Nominal Cartridge Length	Media Formulation
45144	01 = 9 ¾ inch (24.8 cm) Nitrile Gasket 02 = 10 inch (25.4 cm) Nitrile Gasket	1MDS

2. Discs

Basic Disc Number	Disc Size	Media Formulation
64085	01 = 47 mm 02 = 90 mm 03 = 142 mm 04 = 293 mm	1MDS

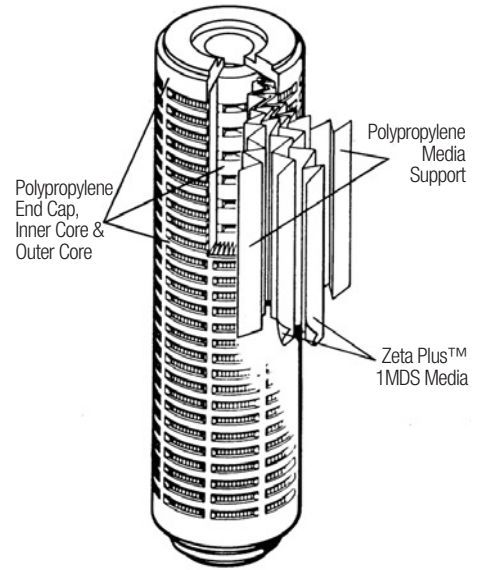

Construction

Zeta Plus 1MDS is an electropositive, charge-modified, glass and cellulose medium. Cartridges are constructed using two layers of pleated medium and are available in 10 inch (25.4 cm) and 9.75 inch (24.8 cm) double open-ended (flat nitrile gasket) configurations. Each cartridge contains approximately 0.39 m² of media and the cartridge may be sterilised using chemical or heat methods (autoclave to 121 °C or 1 hour). When using Zeta Plus 1MDS discs, it is recommended that two layers be used for optimal adsorption. As physical properties of different virus types and environmental conditions vary, it is recommended that the optimum system for the use of 1MDS be established for each laboratory. For more information, contact your local 3M Purification distributor.

Important Notice

3M Purification MAKES NO WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. Since a variety of factors can affect the use and performance of a 3M Purification product in a particular application, some of which are uniquely within the user's knowledge and control, user is responsible for determining whether or not the 3M Purification product is fit for a particular purpose and suitable for user's method of application.

Limited Warranty

3M Purification warrants it this product to be free from defects in material and workmanship during normal use for a period of one (1) year from the date of shipment from the factory. If the Product(s) is (are) defective within this warranty period, your exclusive remedy and 3M Purification's sole obligations shall be, at 3M Purification's option, to replace or repair the Product(s) or refund the original purchase price of the Product(s). This warranty does not apply to failures that result from abuse, misuse, alteration or damage not caused by 3M Purification or failure to properly follow installation and use instructions.

Limitation of Liability: 3M Purification will not be liable for any loss or damage arising from the use of the Product(s), whether direct, indirect, special, incidental, or consequential, regardless of the legal theory asserted, including warranty, contract, negligence or strict liability. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights and you may have other rights which vary from state to state, or country to country.

Zeta Plus and 3M are trademarks of the 3M Company used under license.



LENNTECH

info@lenntech.com Tel. +31-152-610-900
www.lenntech.com Fax. +31-152-616-289

Data may be subject to change without further notice.

v.2 - DOC10180
© 3M 2010. All rights reserved. LITCTZP1MDS1.E - 1110